

Dynamic screen with image and sound
GEWA One – remote control with telephone function

GEWA One is a remote control that helps you manage equipment in your home. You can control TV, music system, receiver for lamps, bed functions, door phone and automatic doors.

In order to keep in touch with your loved ones and quickly reach a support person, you have the opportunity to call and text with GEWA One.

You choose whether you want to use the touch screen, or navigate around with an external control switch. Symbols, colors and voice recordings help you quickly find the right button.

Art.nr 461610



ABILIA

A COMPLETE SOLUTION


GEWA One is a remote control with telephone function, which gives you the opportunity to call and text. You customize how pages and buttons should look and have the opportunity to add speed dials to contacts you want to easily reach.

CONTROL IT YOUR WAY

You only need to tap the screen lightly to activate a button. Select a larger button to get a larger print surface. Instead of having to navigate the remote with an external control switch, Abilia has many different options to choose from.

EASY TO SEE, HEAR AND UNDERSTAND

You can have symbols, different colors and recorded speech on the buttons for easier navigation. Select the number of buttons that a page should display.

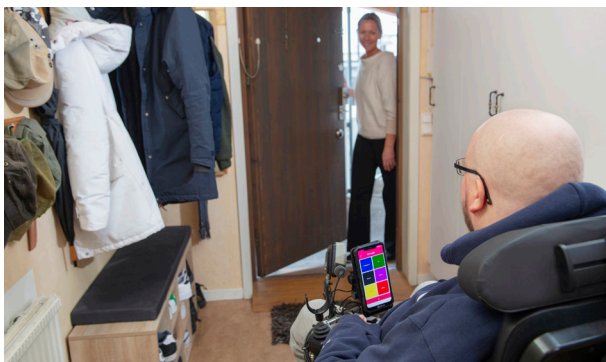

SMART TECHNOLOGY

GEWA One can record signals from eg. TV's usual remote control. GEWA One can also transmit signals to e.g. a GEWA bed receiver for bed control and coded signals for an apartment door. Insert a SIM card for calling and texting with loved ones.

SUPPORT AND SECURITY

When your technician is going to customize the GEWA One to you as a user, this is done smoothly via myAbilia web service. With this service, you as a user have a backup copy of your remote control setup. Your technician can later support you remotely if needed.

EASY TO HANDLE

Accessories such as a holder specially designed for the GEWA One allow you to see and manage the remote control from both wheelchair and bed.

DELIVERED WITH

- Mobile unit, restricted to GEWA One
- Charger mobile unit incl USB-cable
- GEWA Multibox
- Charger GEWA Multibox incl USB-cable
- Headphones, 3,5 mm incl USB-C adapter
- USB-cabel USB-C to USB-A
- SIM-tray installation tool
- User manual
- Voucher with license key

TECHNICAL INFORMATION GEWA ONE

MOBILE UNIT

Screen size	6.5"
Size	178 x 84.7 x 12,9 mm
Weight	260 g
OS	Android 10

MULTIBOX

Size	108 x 67 x 22 mm
Weight	153 g
Connections:	3.5 mm, Mikro-USB
Control signals:	IR, Radio

ACCESSORIES

Holder	461186
Miniarm	451030
Universalarm	621966

ART NR

