



For a social life
Lightwriter® series

The Lightwriter SL40 is a dedicated device for those who are literate, but unable to communicate effectively through speech. With a Lightwriter you can participate in social interactions while at work or during leisure time, even in noisy, sunny outdoor or poorly lit environments.

Using dual screens, your communication partner can read the text whilst maintaining a natural face-to-face dialogue, including engagement through gestures and body language.



ABILIA
Toby Churchill




The Lightwriter is a communication tool where you just press the On button, start typing, and off you go. The word prediction enables you to type more effectively. The Lightwriter can be used both for face-to-face communication and remote communication via SMS and mobile phone.

FOR A MORE ACTIVE LIFE

Do you want to socialise with your friends and be able to communicate with them in different environments and settings? The Lightwriter gives you access to technology specifically designed for communication, that can help to enrich your life. It will encourage and support engagement with your family, friends and the community.

STAY IN TOUCH WITH YOUR FRIENDS AND FAMILY

To stay in touch with friends and family is important for everyone. With a SIM card inserted in your Lightwriter SL40 you can send text messages, plus the SL40 Connect also allows you to make mobile phone calls. Ideal for keeping friends and family updated when you're out for the day.

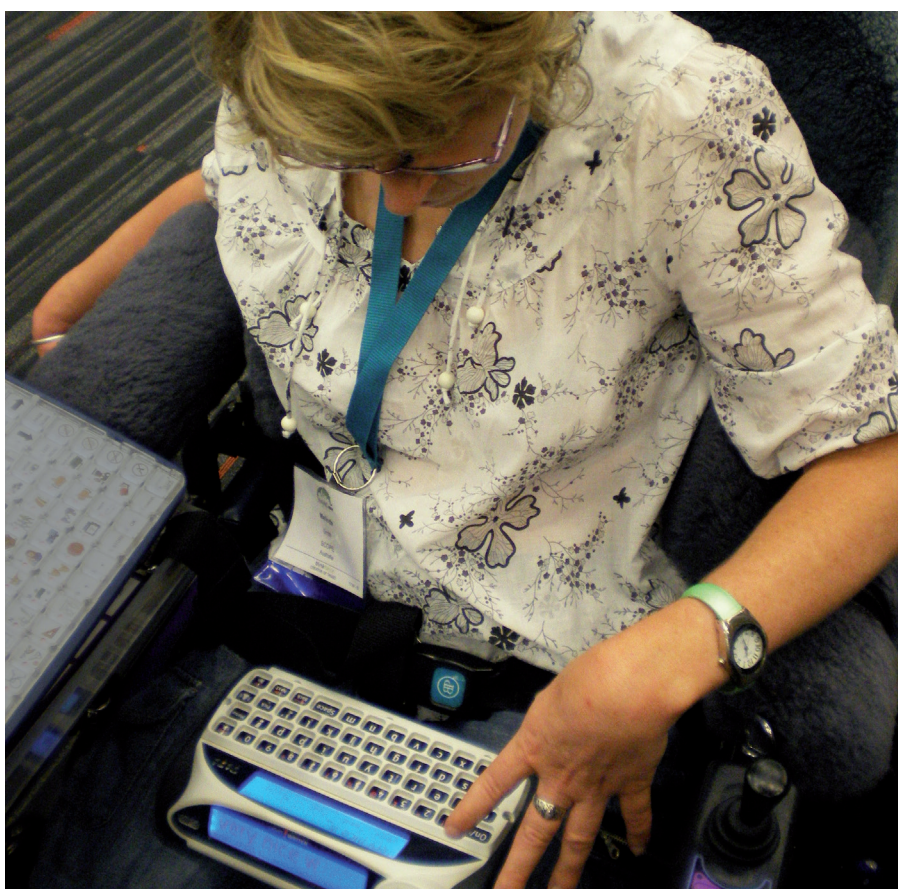

EASY COMMUNICATION

The Lightwriter has dual screens, allowing your communication partner to read the text whilst maintaining a natural face-to-face dialogue, including engagement through gestures and body language. The built-in speaker has very good sound quality and faces your communication partner.

The Lightwriter can predict which words you are going to use and proposes various alternatives. It is possible to pre-store words and sentences for even faster communication. With the SL40 Connect you can back-up your personalised vocabulary to a USB stick or a PC.



"I like that I can talk to many different people when I'm out"

TAKE IT EVERYWHERE

The Lightwriter is small, portable and robust. You can take it everywhere as it is designed to withstand the bumps and knocks of everyday life. Long battery life means you can use your Lightwriter as much as you want all day every day, if you recharge it every night.

EASY TO USE

The Lightwriter is easy to use – just press a button and you are ready to go! It has a user friendly menu and excellent sound quality using text-to-speech technology (Acapela).

Key sensitivity can be adjusted and each Lightwriter comes with a set of keyguards. These help to guide your fingers to the correct location and avoid accidental key presses.

FLEXIBLE ACCESS OPTIONS

The Lightwriter SL40 is also available with scanning. Scanning is useful for those who do not have the physical ability or energy to use a full keyboard. The keys will light up row-by-row or key-by-key and are selected by activating a switch.

Switches can be activated by a variety of methods including pressing a switch by hand, foot or head, or operating a sip/puff switch by mouth. They can also differ in pressure sensitivity, size and texture in order to suit an individual's needs.



WANT TO KNOW MORE?

Additional product features can be found in the user guide on our website.

ABOUT ABILIA


Abilia gives people the power to become independent and take control of their daily lives. Based on more than 40 years of experience and supported by research evidence, we know that our assistive technology helps people to have structure over their day, communicate, control appliances in their home and call for attention. We believe that everyone should have the same opportunities to participate in society on their own terms, and are passionate about enabling people access to a more fulfilling life.

CONTACT US

Phone: 01954 281 210

Email: info@abilia.uk

Visit www.toby-churchill.com for more information.

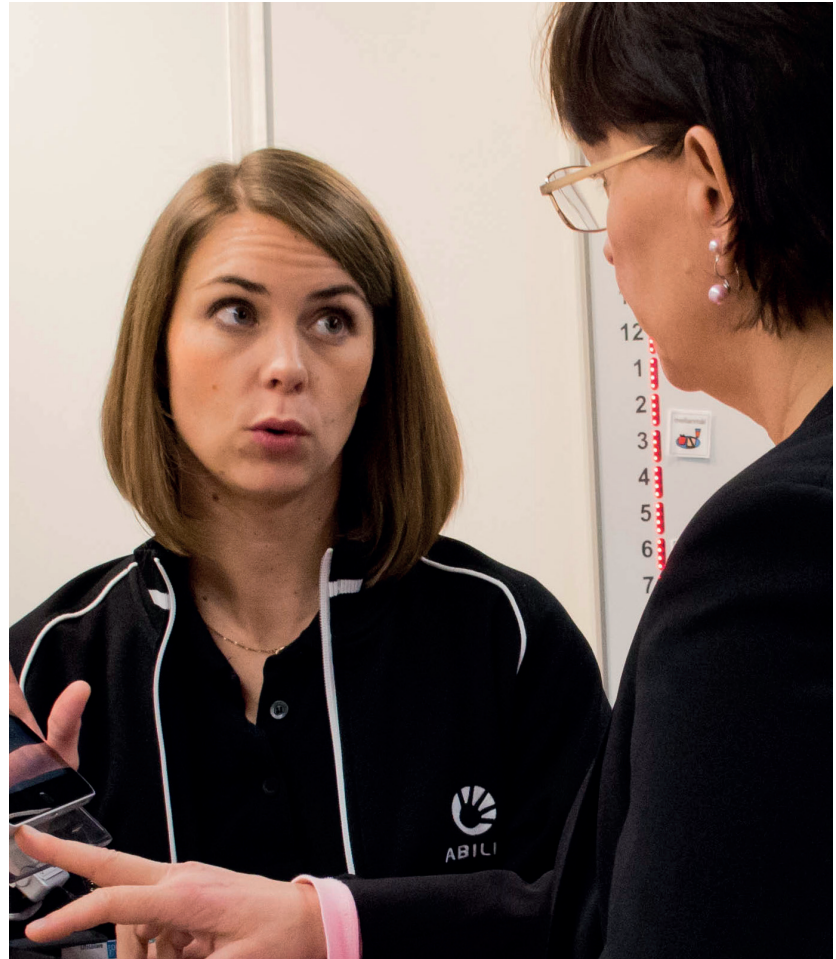

PRODUCT MODELS

Lightwriter SL40.

Lightwriter SL40 Connect.

Lightwriter SL40 Scanner.

LIGHTWRITER COMES COMPLETE WITH

Three different keyguards, carry bag, SIM Card, charger, tool kit, user guide and getting started leaflet, memory stick (Connect model only).

ACCESSORIES

In Use Bag.

Trabasack Mini - Slim Lap Tray Bag.

Complete Hand, Neck & Waist Strap Set.

In-car charger.

TECHNICAL SPECIFICATIONS

Size: 240mm x 140mm x 55mm

Weight: 910 g

Speech synthesis: Acapela HQ.

Battery: Rechargeable nickel metal hydride.

Battery life: All day battery life with standard usage.